


D'z IMAGE Viewer C

Setting Up the Wi-Fi Connection on your PC

- From "Start" **1** on the Windows desktop screen, click "Settings" **2**.



- Click "Network & Internet" **3**.

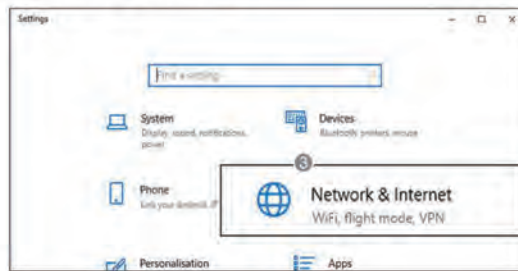

- Click "Wi-Fi" **4** and check that the Wi-Fi setting is enabled.
- Click "Show available networks" **5**.



- Select SSID **6**.

- Wi-Fi connection via AP mode

Select the camera SSID from the PC network list (SSID) and click "Connect".

SSID : DZ-C100
Password: 00000000

- ※ SSID and password are default settings.
- ※ To protect personal information, it is recommended that you change the password.

- Wi-Fi connection via STATION mode

Select the router SSID from the PC network list (SSID) and click "Connect".

