



D'z IMAGE Viewer D

Setting Up the Wi-Fi Connection on your PC


- From “Start ①” on the Windows desktop screen, click “Settings ②”.

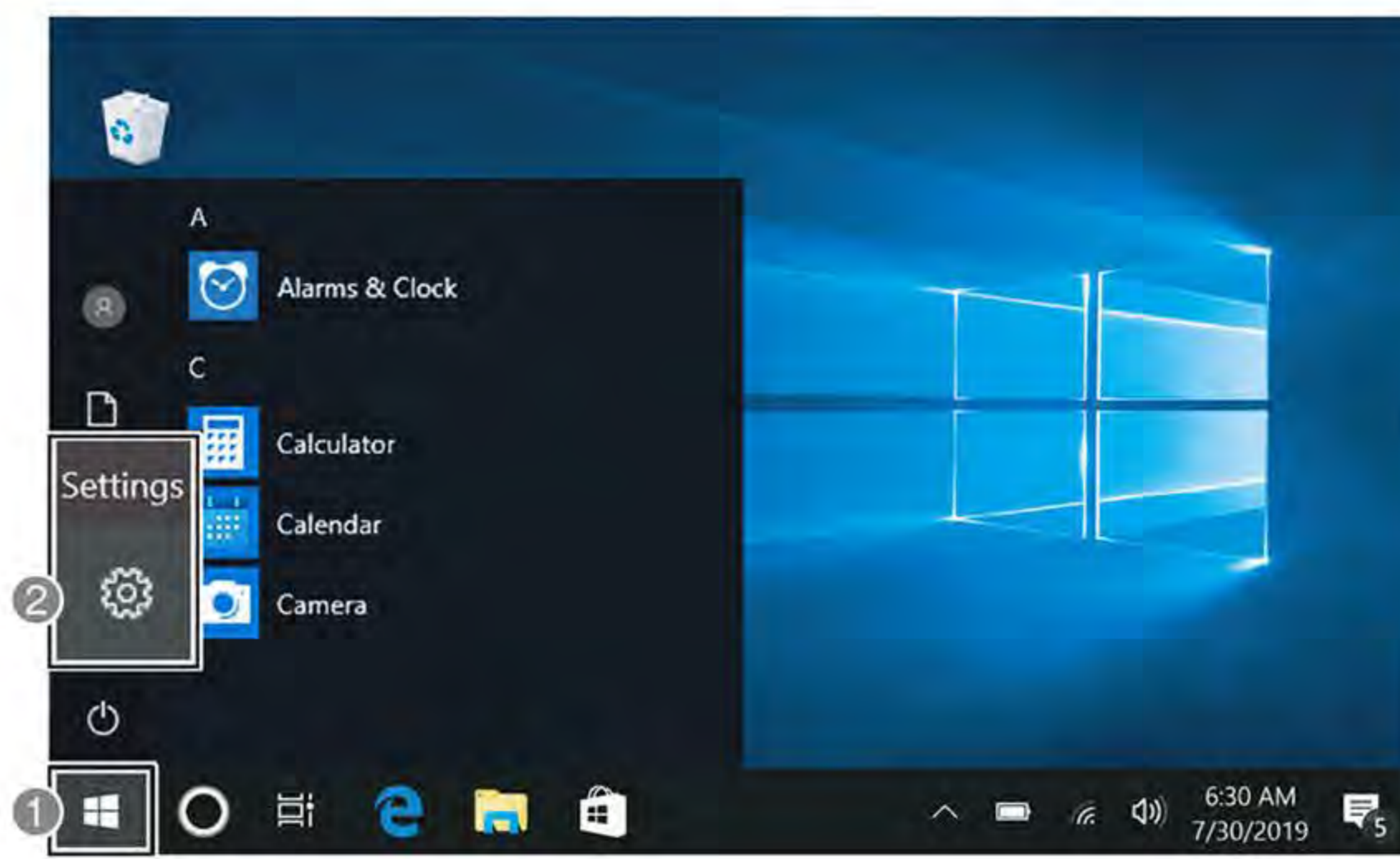

- Click “Network & Internet ③”.

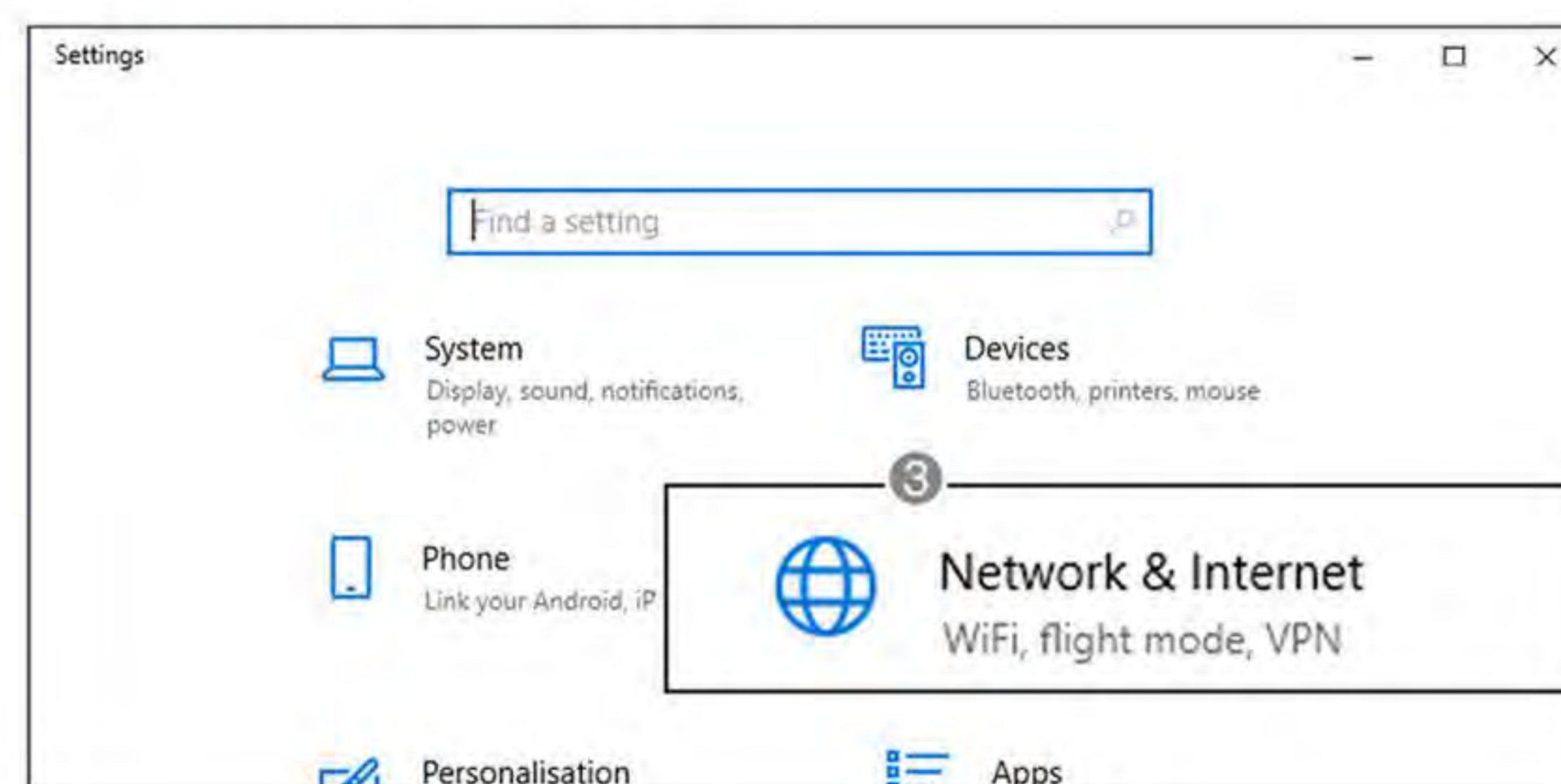

- Click “Wi-Fi ④” and check that the Wi-Fi setting is enabled.
- Click “Show available networks ⑤”.

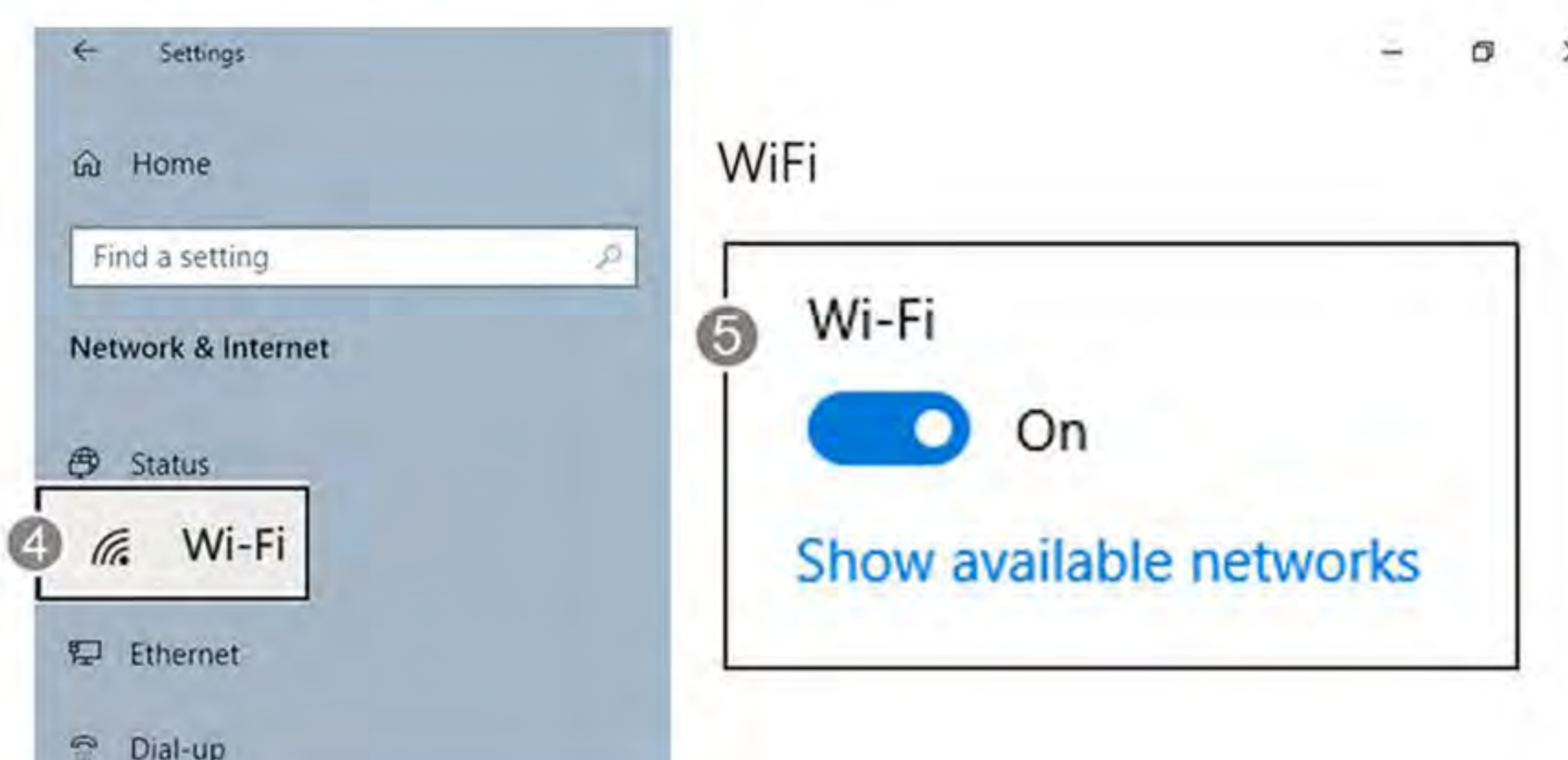

- Select SSID ⑥.

Wi-Fi connection via AP mode

Select the camera SSID from the PC network list (SSID) and click "Connect".

SSID : DZ-D100
Password: 00000000

- ※ SSID and password are default settings.
- ※ To protect personal information, it is recommended that you change the password.

Wi-Fi connection via STATION mode

Select the router SSID from the PC network list (SSID) and click "Connect".

